

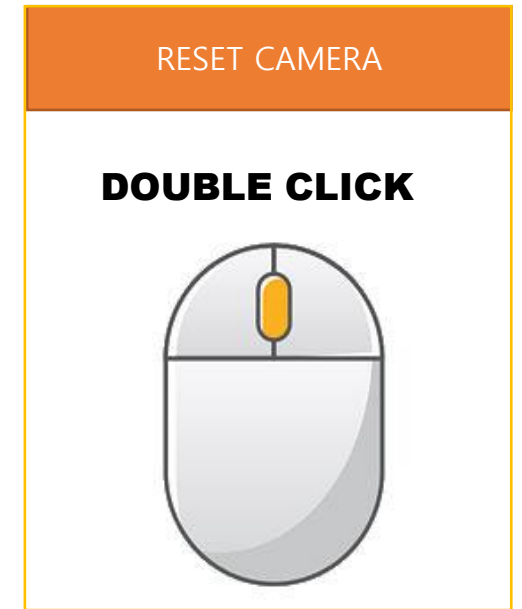
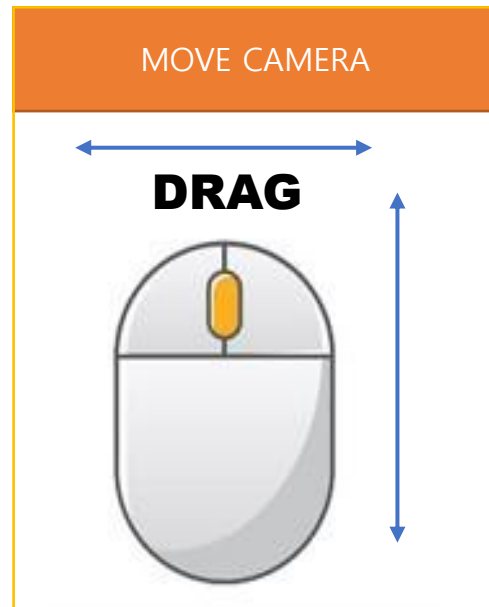


# COMMON CAMERA CONTROL

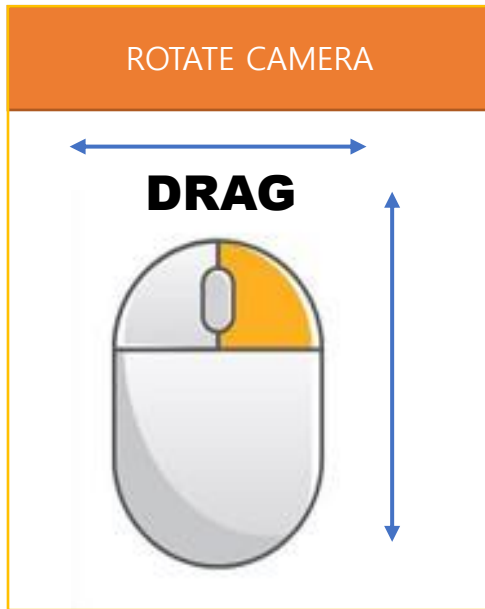

# PROTOCOL EDITOR



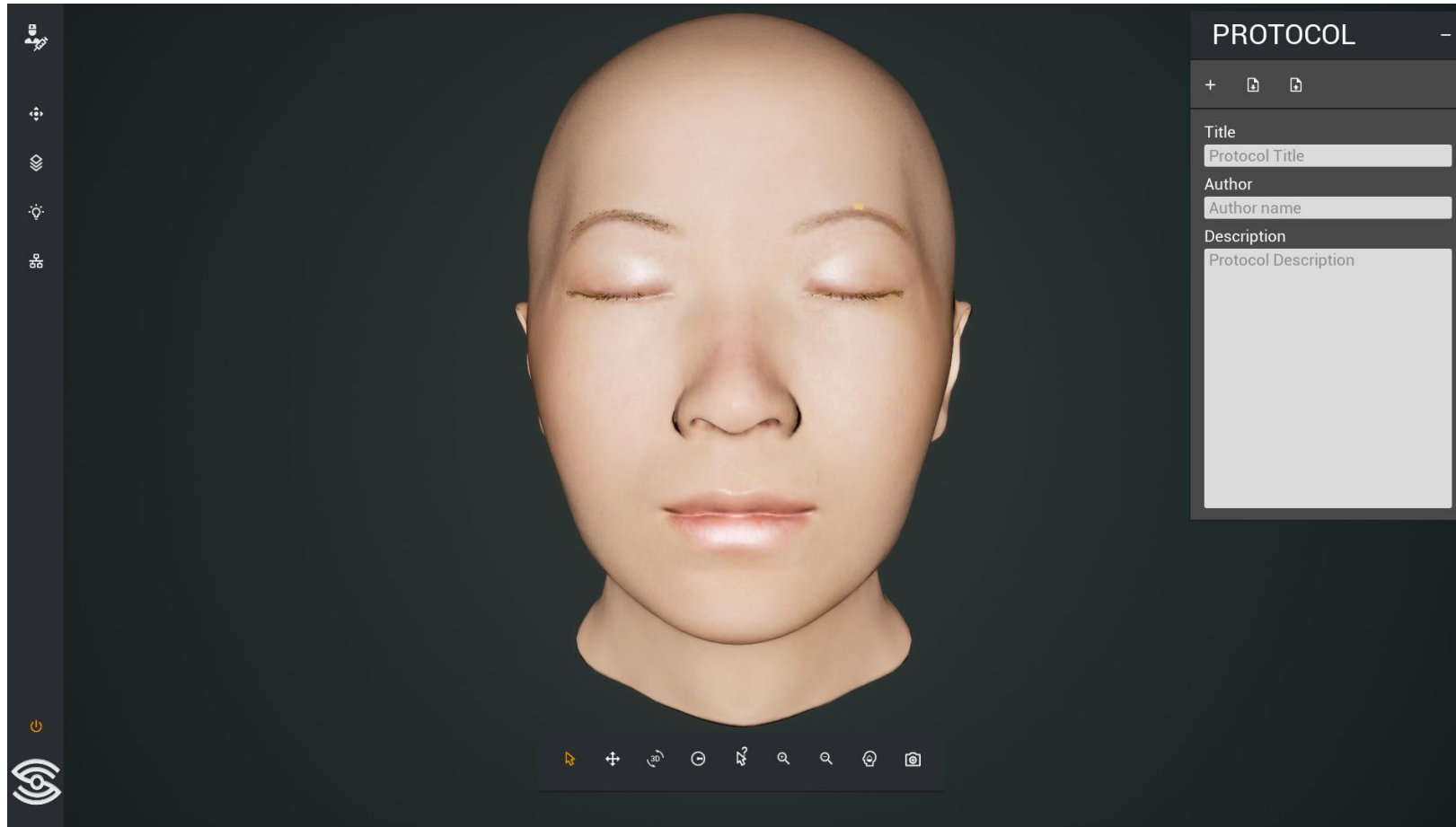


# Side Menu Bar


- 0. Start Training
- 1. Camera Menu
- 2. Layer View Menu
- 3. Highlight View Menu
- 4. Anatomy Tree Menu




# Camera Menu

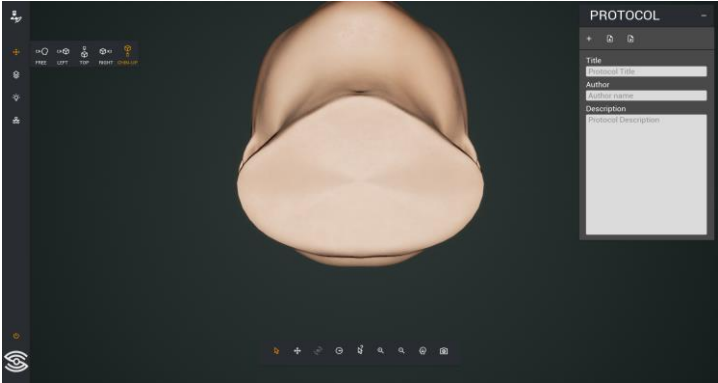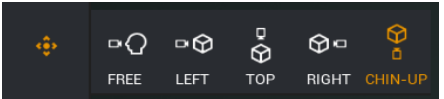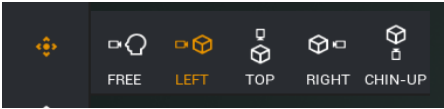




# Layer View Menu

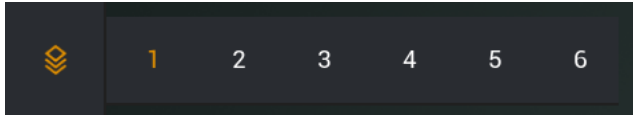


1. Skin : All Anatomy
2. Superficial Fat Layer
  - Muscle, Superficial Fat, Vessels, SMAS, Catilage, Skull
3. Deep Fat Layer
  - Deep Fat, Vessels, Gland, Skull
4. Muscle Only
  - Muscle, Vessels, SMAS, Catilage, Gland, Skull
5. Vessels
  - Vessels, SMAS, Catilage, Gland, Skull
6. Nerve Only
  - Nervous, Gland, Skull

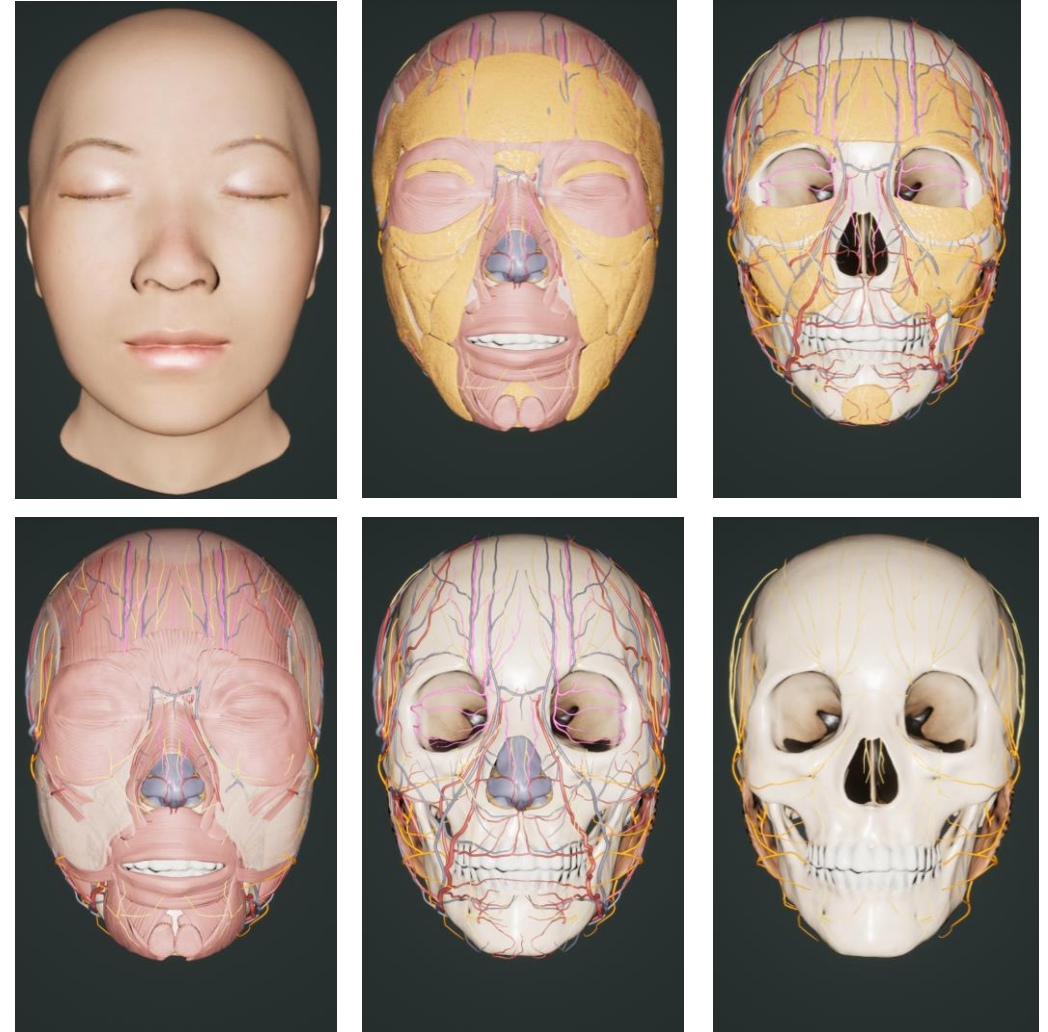

# Highlight View Menu

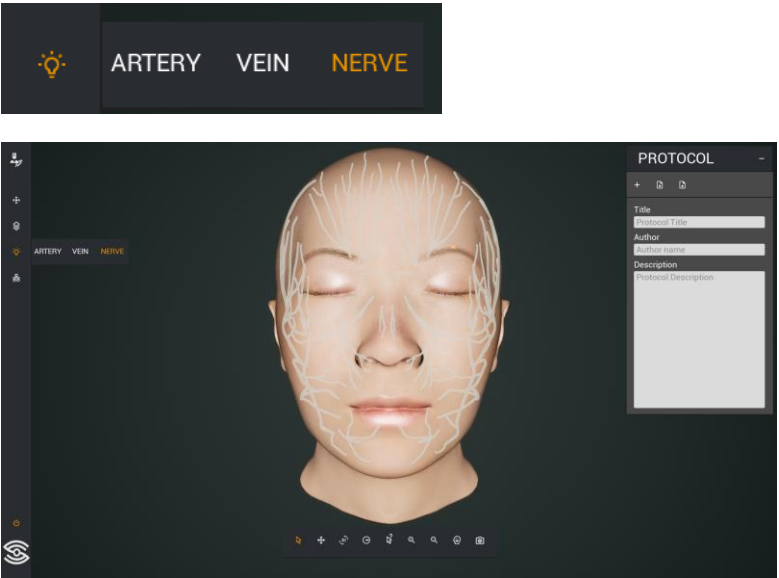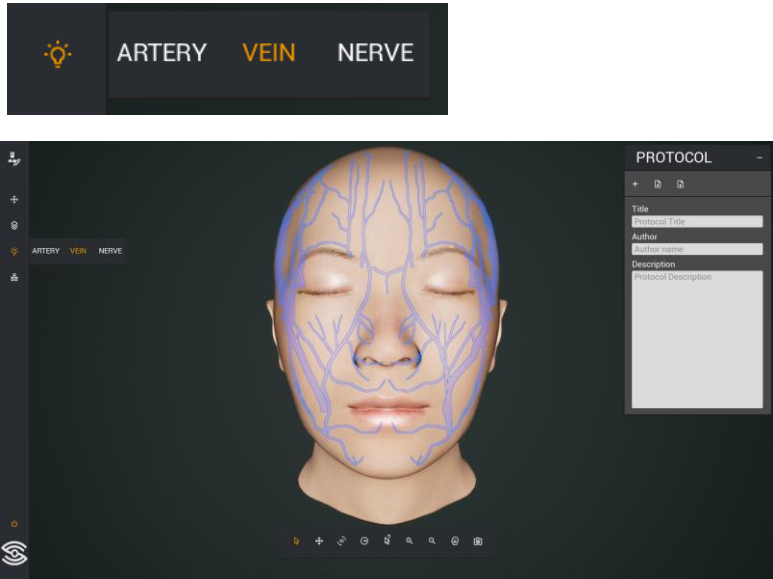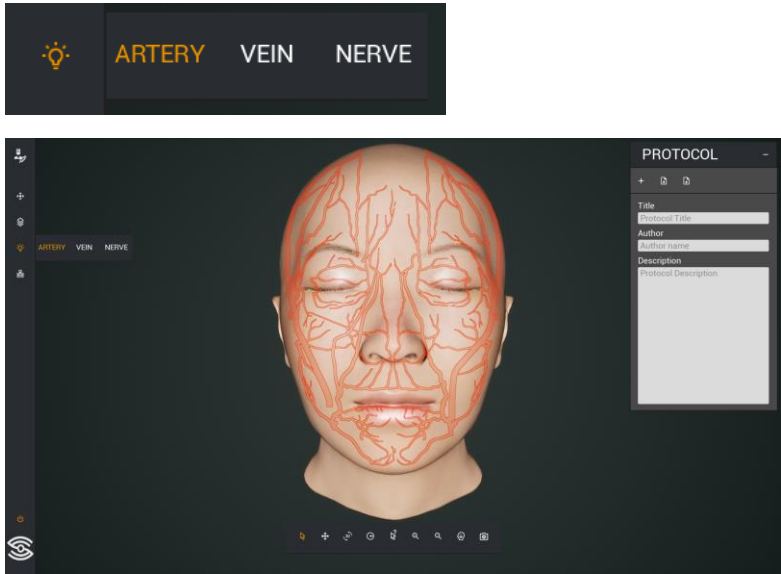


# Anatomy Tree

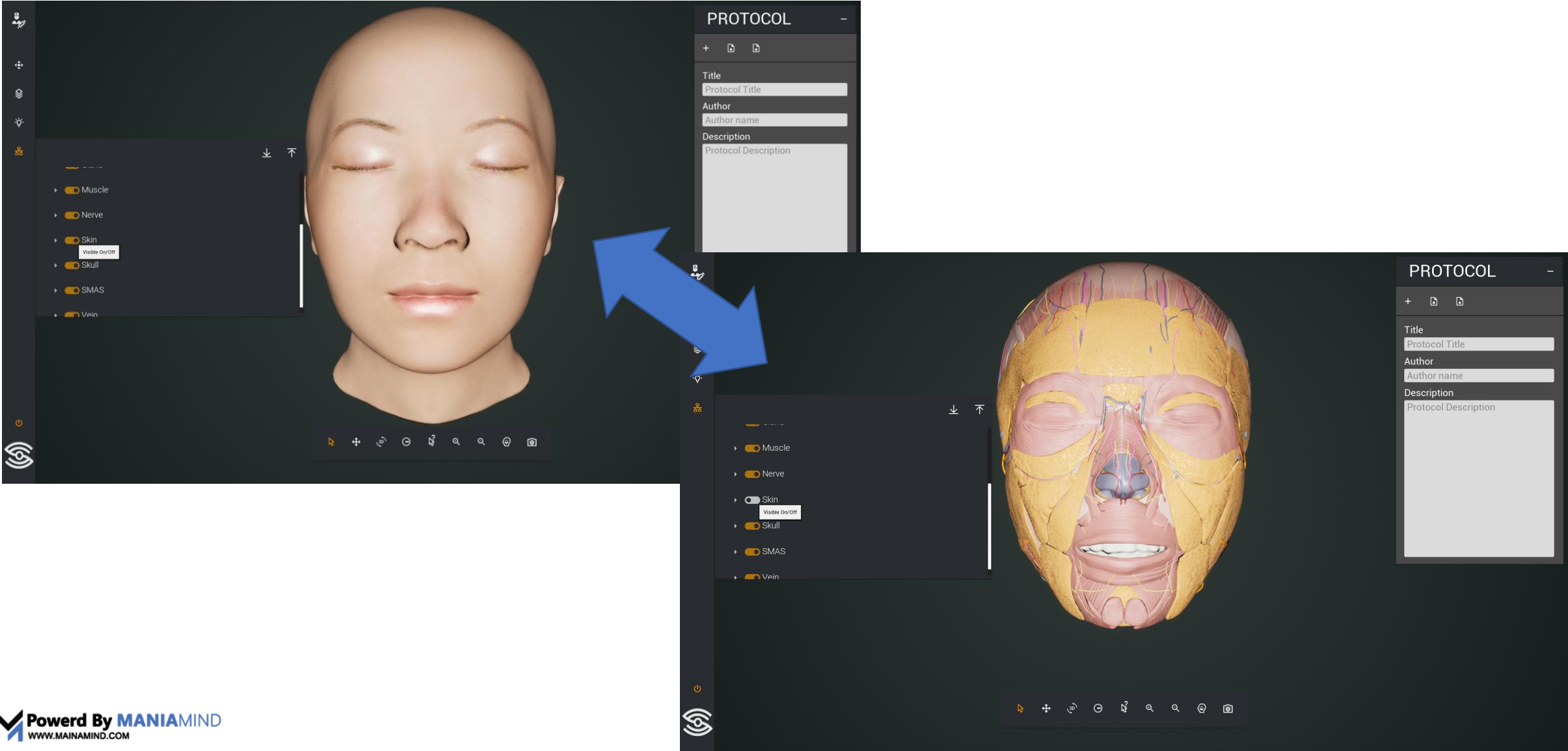




# Bottom Tool Bar

1. Select Tool
2. Move Camera
3. Rotate Camera
4. Create Protocol Point
5. Toggle Anatomy Label
6. Zoom-In
7. Zoom-Out
8. Toggle Skin Translucent Mode
9. Screenshot

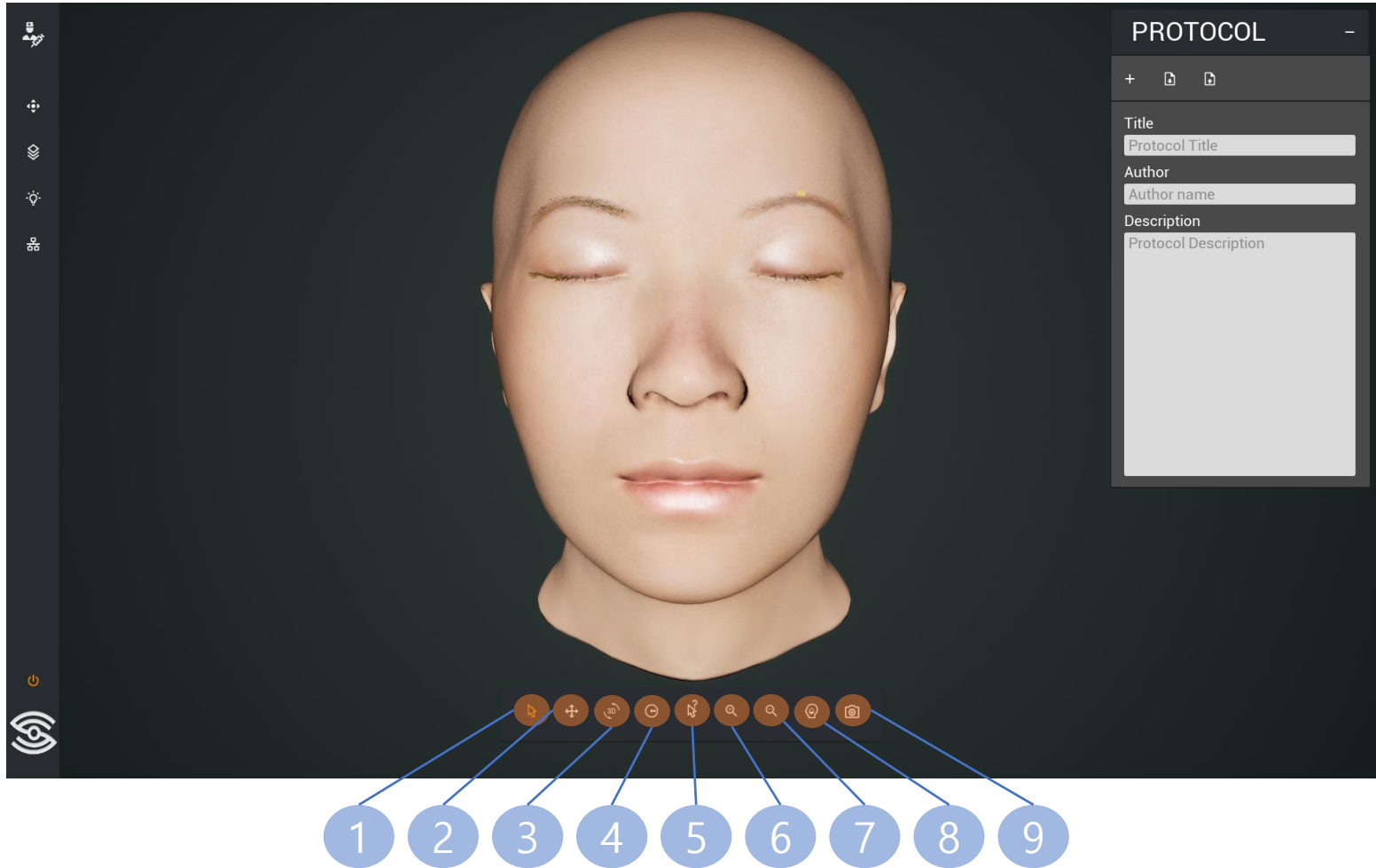




# Move/Rotate Camera



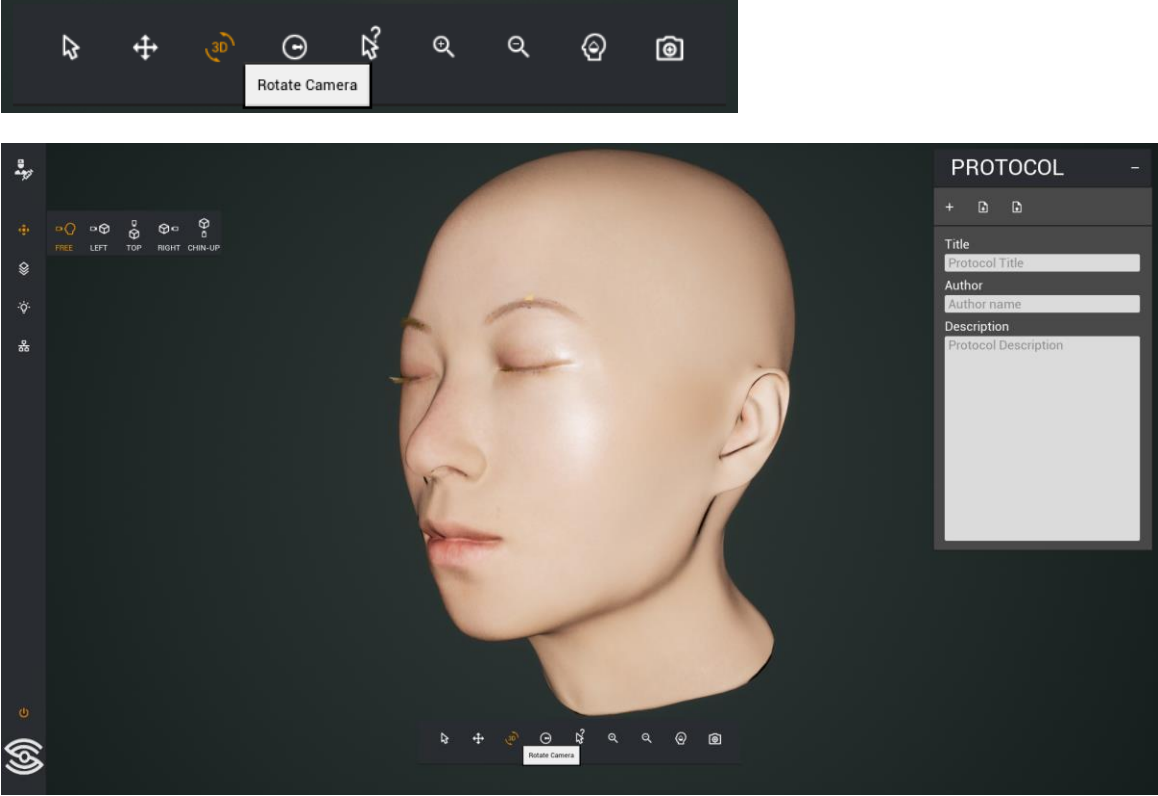


# Toggle Anatomy Label


- Show anatomy name tooltip on mouse pointer.
- It works on "Select" mode.

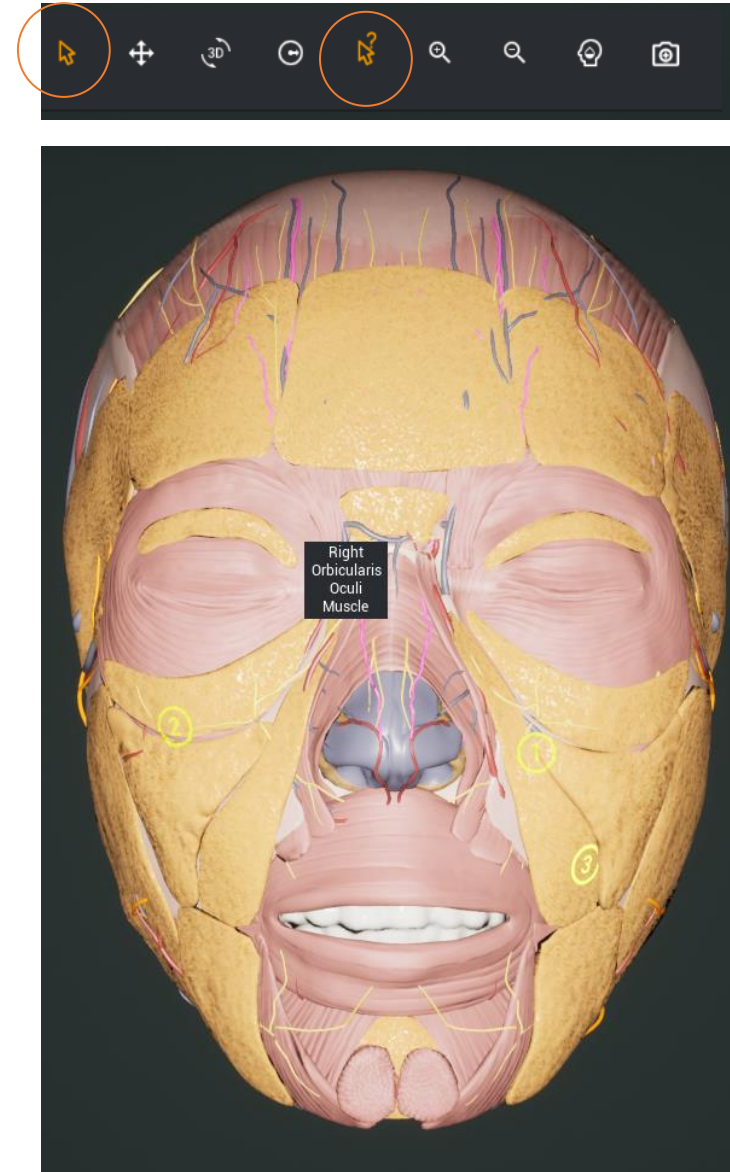

# Toggle Skin Translucent

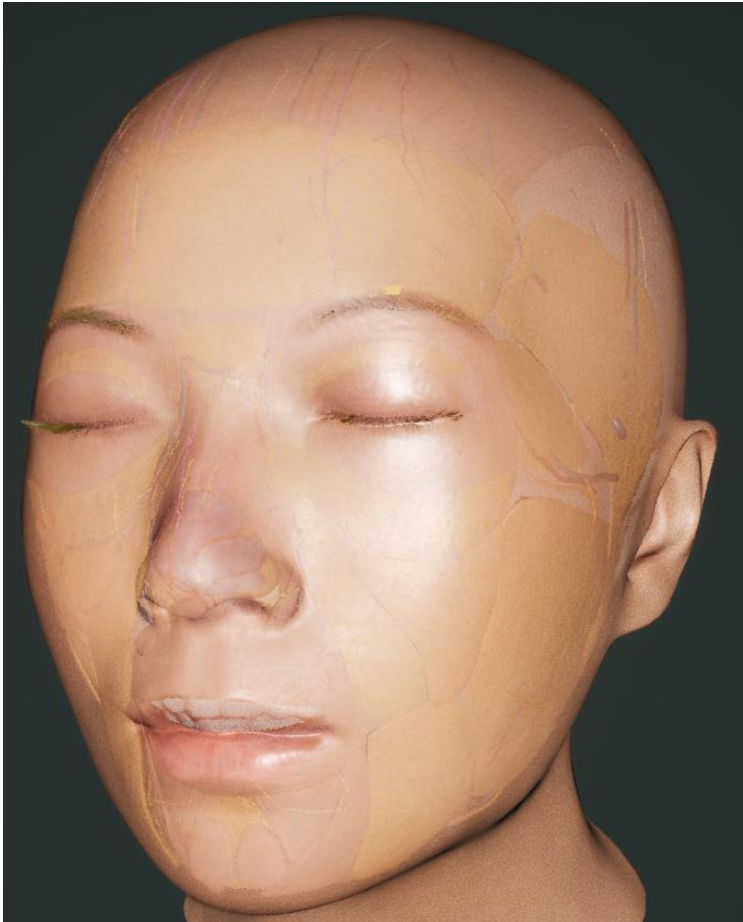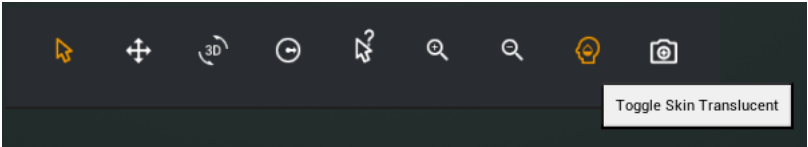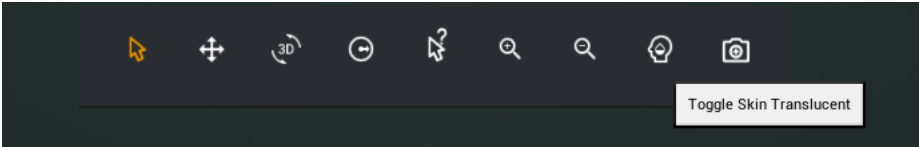


# Screenshot

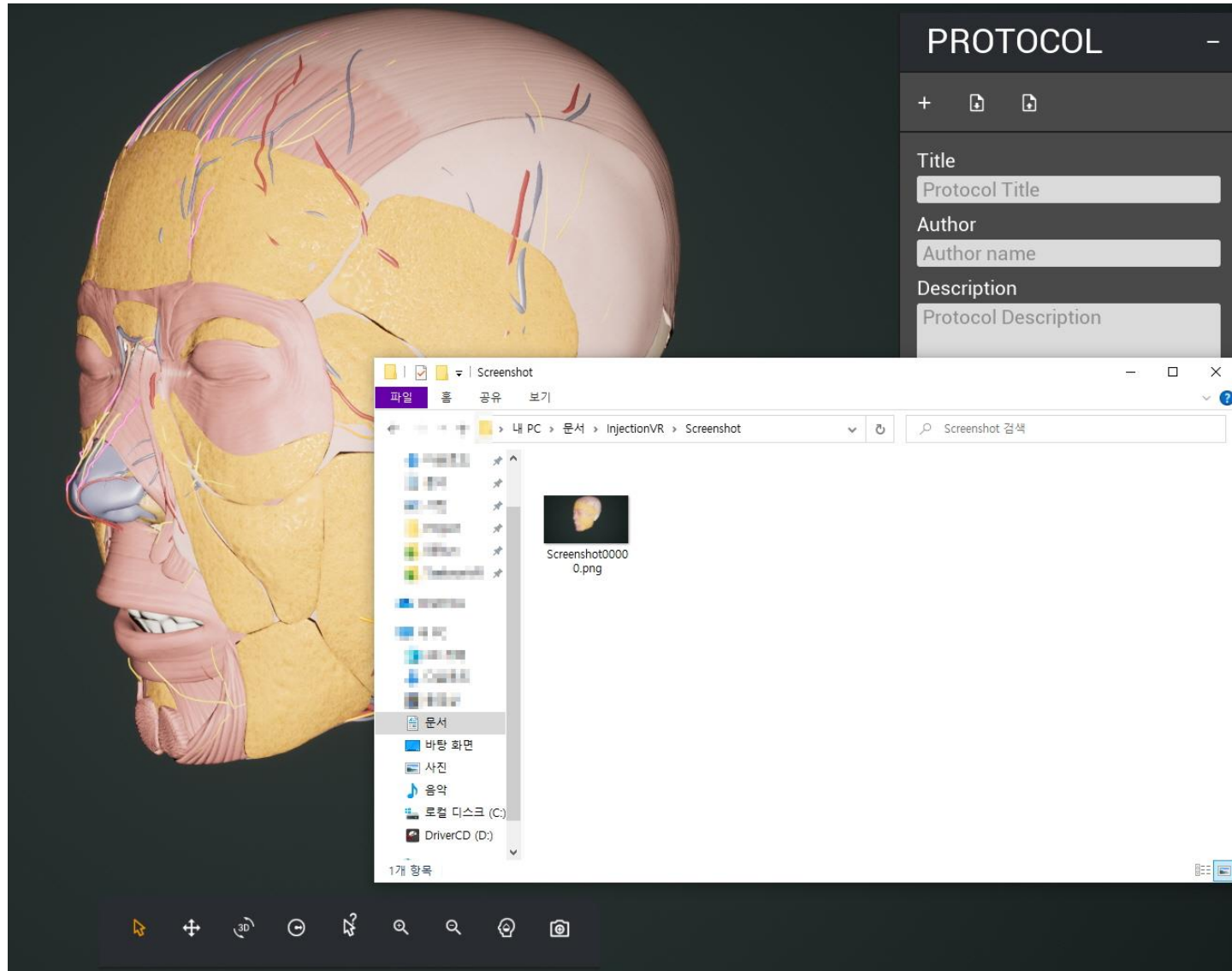

Screenshot path : My Document/InjectionVR

# Protocol Editor Panel


1. New protocol button
2. Load protocol button
3. Save current protocol button
4. Protocol title
5. Protocol author
6. Protocol description

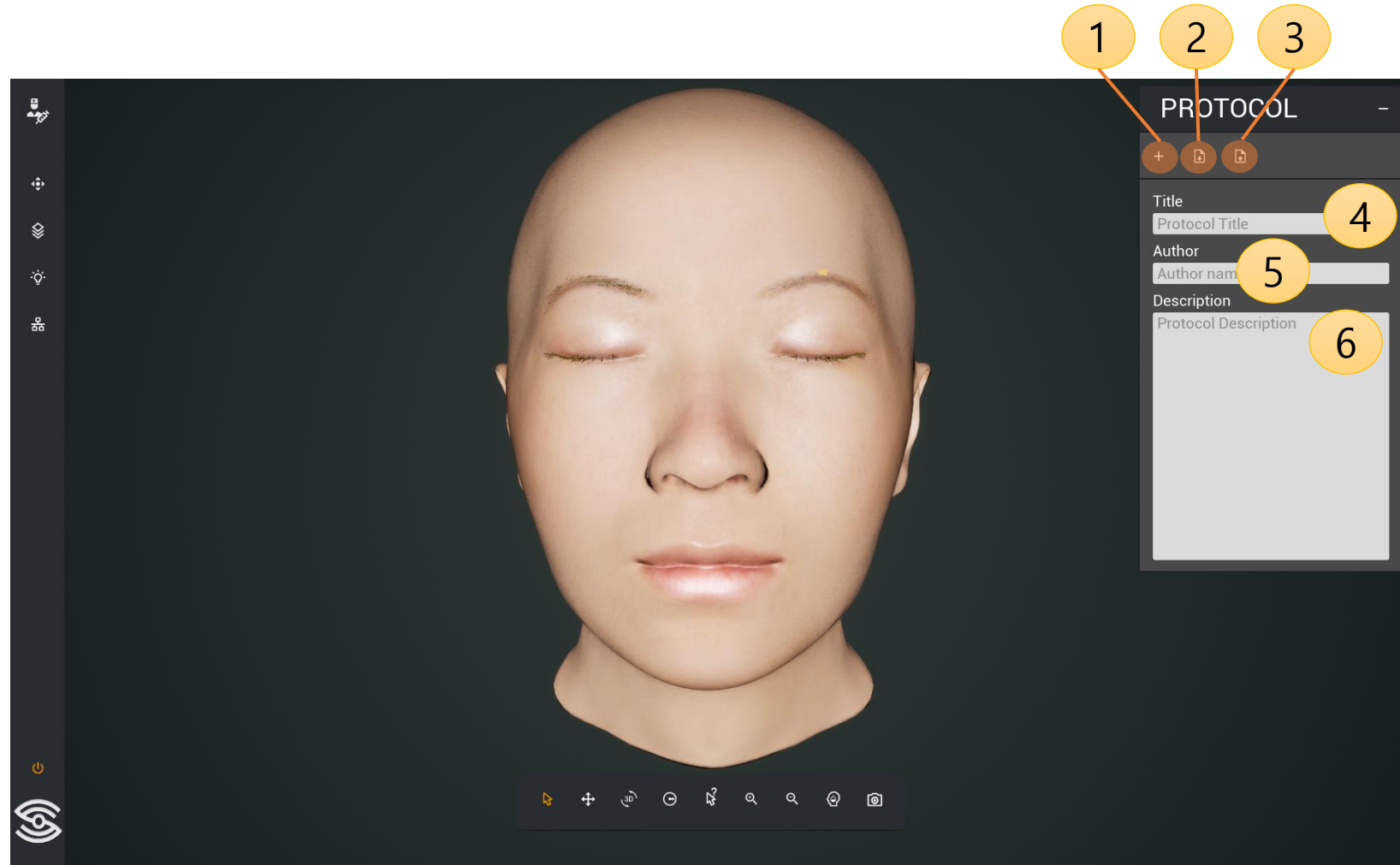

# How to create new protocol

1. Click [New Protocol] button
2. Fill [Protocol Title]
3. Input author name on [Protocol Author]
4. Input description on [Protocol Description]
5. Create / Move points (Details in next pages)
6. Click [Save Protocol] button

Protocol file save path : InstallPath/InjectionVR/Protocol/  
Default install path : C:/SurgicalMind/InjectionVR

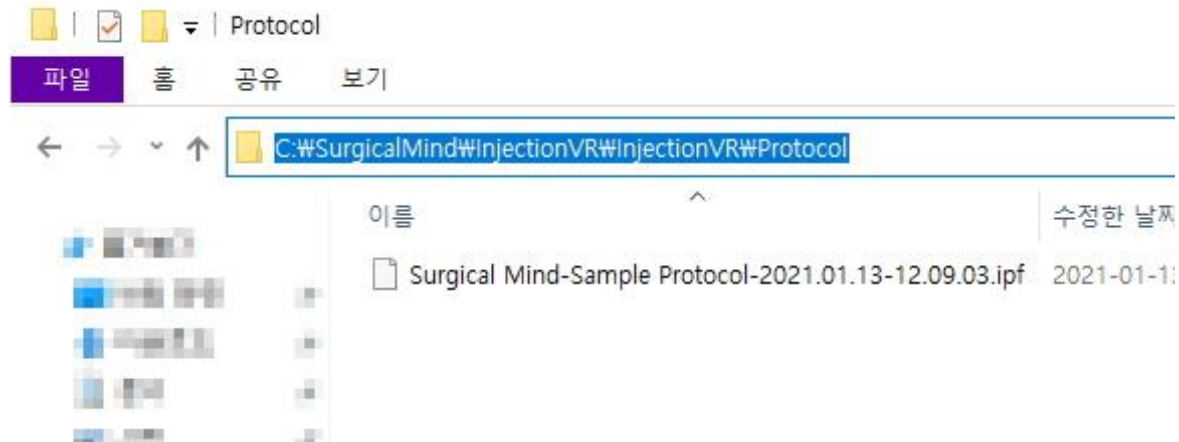

## PROTOCOL

+

Title

Sample Protocol

Author

Surgical Mind



Description

Sample Protocol



# Create Protocol Point

1. Select "Create New Point" Tool on the Bottom Tool Bar
2. Click left mouse button anywhere on the skin.

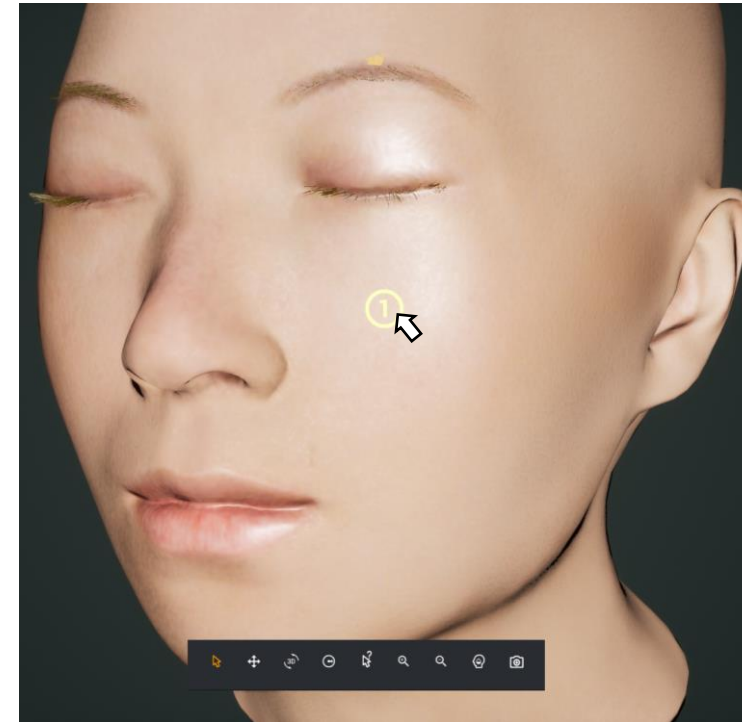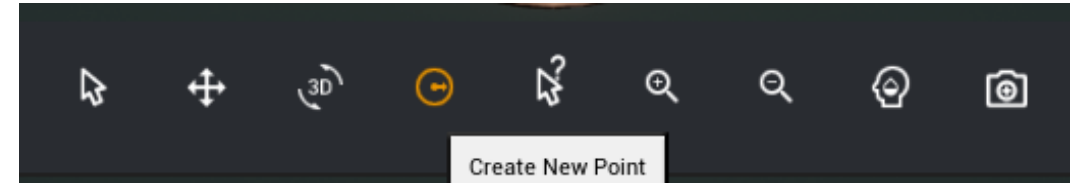



**CLICK**





# Select & Move Point


1. Select "Select" tool on the bottom tool bar
2. Select any point
3. Press left mouse button and drag

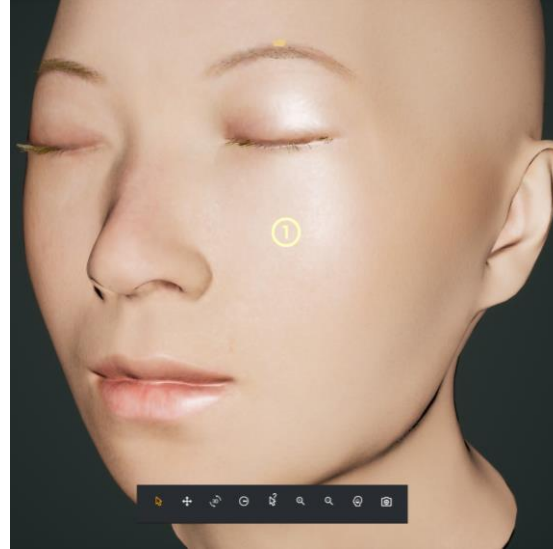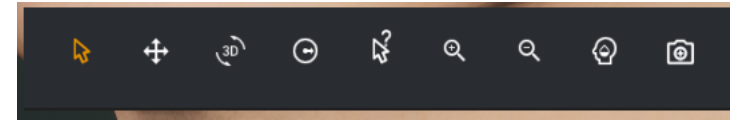

# Delete Protocol Point


1. Select "Select" tool on the bottom tool bar
2. Move mouse and select a point
3. Click right mouse button
4. Click trash icon on the popup

Also you can set the product amount at point.

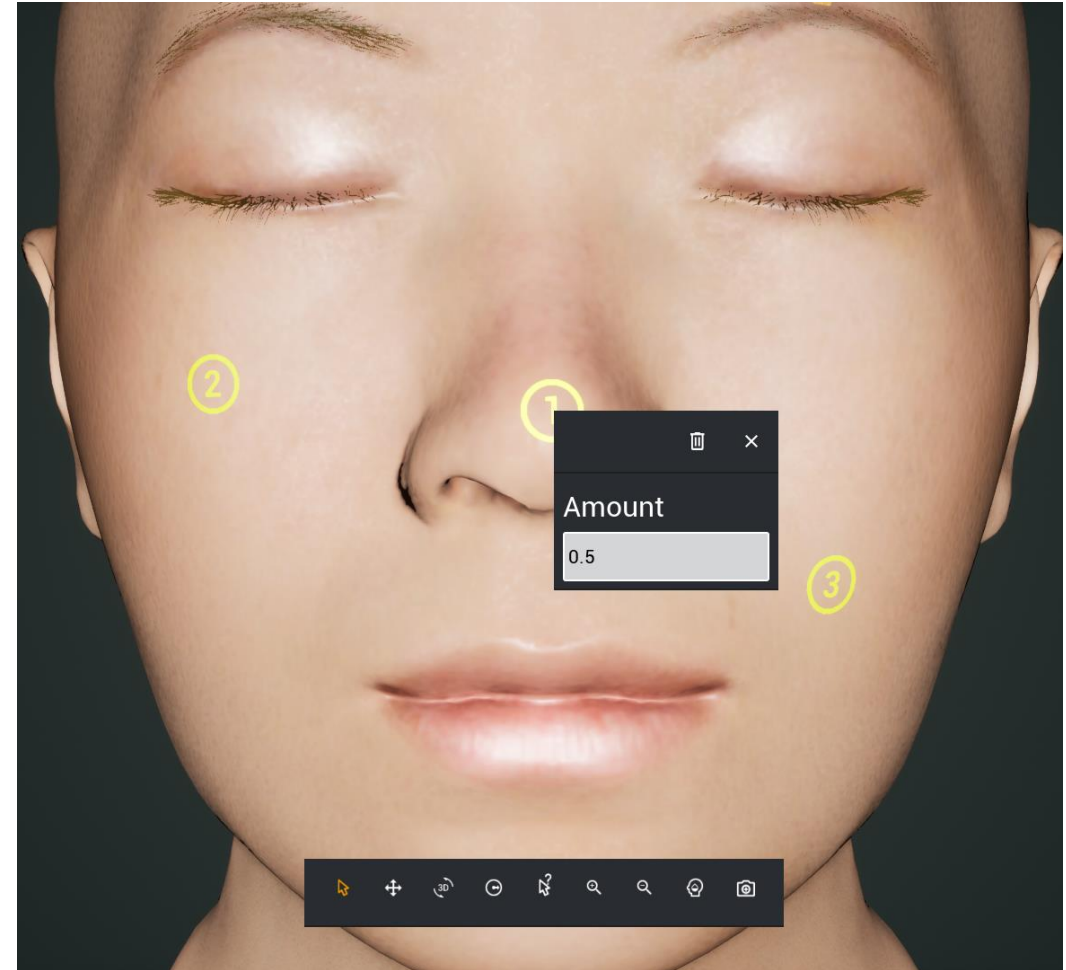

# Load Protocol from file

1. Click [Load Protocol] button on the Protocol editor panel
2. Select protocol to load
3. Click [Load] button

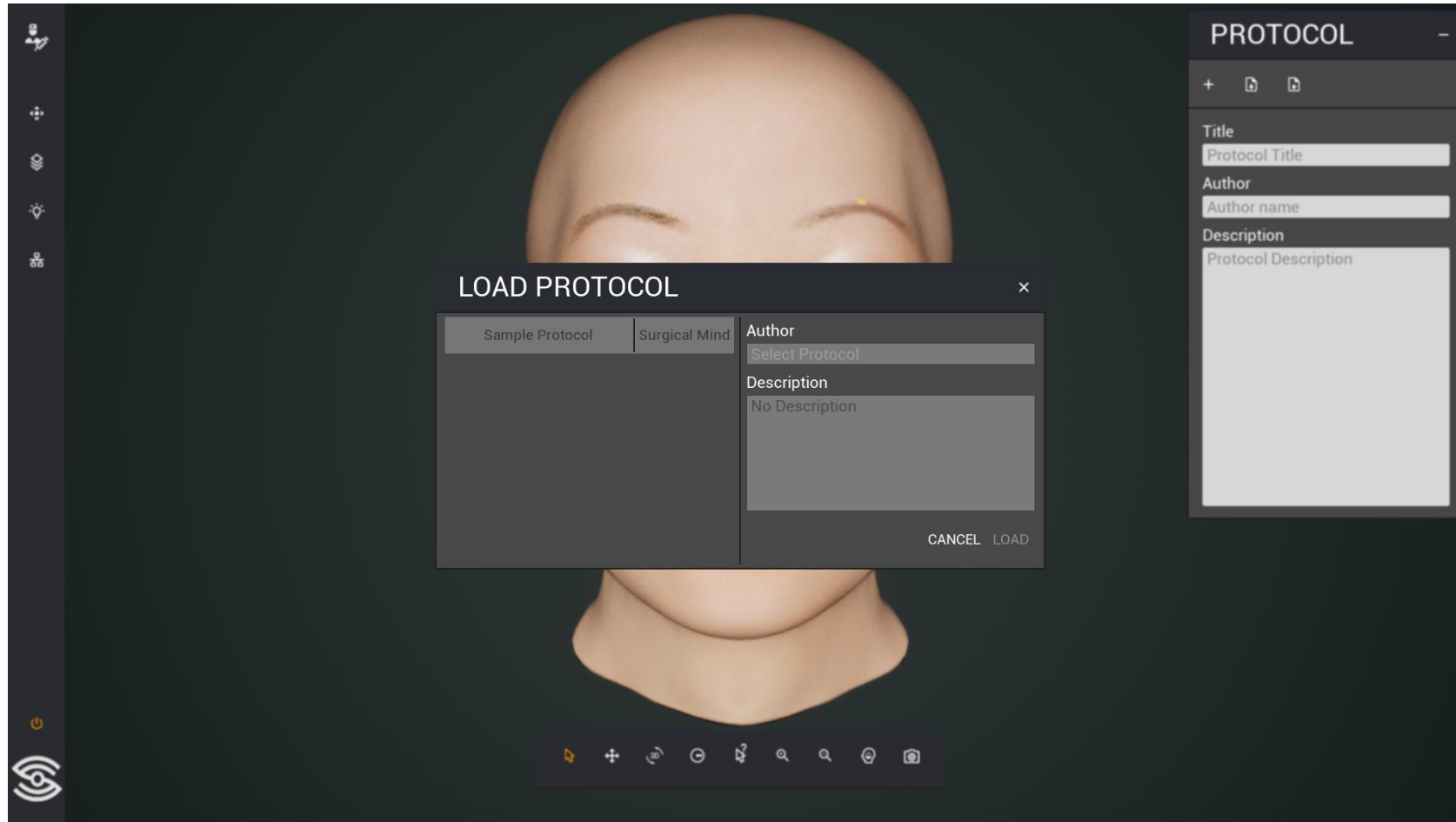


# TRAINING




# Toggle Menu Bar

Mouse click on table to toggle menu bar

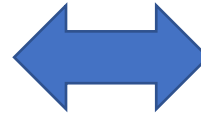

# Training Menu Bar

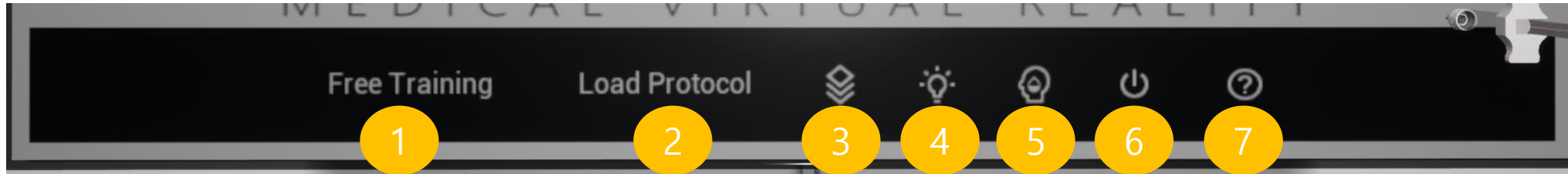

1. Free Training (Injection training without protocol point marker)
2. Load Protocol
3. Layer View Menu
4. Highlight View Menu
5. Toggle Skin Translucent
6. End Training (Go to Protocol Editor)
7. Show Information pop up